

A'TOMIC PULSE™

VERSATILE MINI-FRAGMENT FIXATION



**SMOOTH
PROFILE**

**ASTONISHING
STRENGTH**

Why
put a
Square Peg
in a
Round Hole?

RMA
ORTHO

Implant shown: **A'TOMIC PULSE™** 13X10mm



Zero-Height Release.

Designed to prevent disruption of meticulous reductions.

Ultra-Low Profile.

When troughed, our unique bridge cross-section (blue) allows for sub-cortical placement while maximizing the strength and stiffness of the implant bridge.



Choose implant depth:
1) Tamp only, 2) 1mm trough, 3) 2mm trough


Precision Instruments.

Cannulated drilling option for precise implant placement.




Strength Where it Counts.

Patent-pending design with round legs maximizes implant strength, stiffness, and dynamic compression while minimizing strains that can compromise implant integrity.

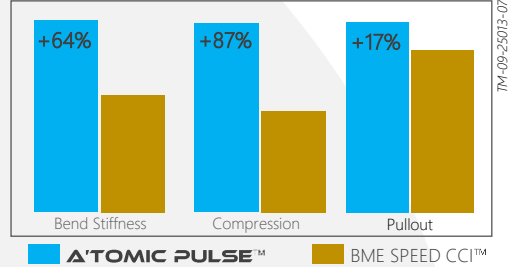

Silky Smooth.

Smooth contours designed to protect against soft tissue irritation caused by hardware.

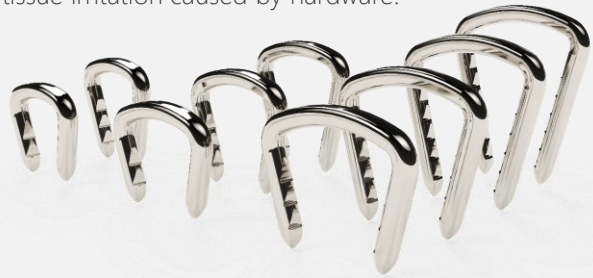

SKU	DESCRIPTION
A'TOMIC PULSE™	
07.20.1108.21A	A'TOMIC PULSE 2.0MM NITI IMPLANT KIT 11X08MM 2-LEG, STERILE
07.20.1110.21A	A'TOMIC PULSE 2.0MM NITI IMPLANT KIT 11X10MM 2-LEG, STERILE
07.20.1308.21A	A'TOMIC PULSE 2.0MM NITI IMPLANT KIT 13X08MM 2-LEG, STERILE
07.20.1310.21A	A'TOMIC PULSE 2.0MM NITI IMPLANT KIT 13X10MM 2-LEG, STERILE
07.20.1312.21A	A'TOMIC PULSE 2.0MM NITI IMPLANT KIT 13X12MM 2-LEG, STERILE
07.20.1510.21A	A'TOMIC PULSE 2.0MM NITI IMPLANT KIT 15X10MM 2-LEG, STERILE
07.20.1512.21A	A'TOMIC PULSE 2.0MM NITI IMPLANT KIT 15X12MM 2-LEG, STERILE
07.20.1514.21A	A'TOMIC PULSE 2.0MM NITI IMPLANT KIT 15X14MM 2-LEG, STERILE
07.20.1516.21A	A'TOMIC PULSE 2.0MM NITI IMPLANT KIT 15X16MM 2-LEG, STERILE
ANCILLARY INSTRUMENTS	
07.20.IPK	A'TOMIC 2.0MM INSTRUMENT PROCEDURE KIT, STERILE
07.20.CDK	A'TOMIC 2.0MM CANNULATED DRILL KIT, STERILE
07.20.QDS	A'TOMIC 2.0MM QUICK DISCONNECT DRILL, STERILE (SINGLE-USE)



RMR Ortho, LLC
4600 Lockhill-Selma Rd,
Suite 107
San Antonio, TX 78249
Phone: (888) 476-7210
Email: customersupport@rmrortho.com
www.rmrortho.com