



G&G TECHNOLOGIES, INC.

G&G TECHNOLOGIES, INC. INTRODUCES THE

G2 VERSAFLOW PUMP


FOR HAZARDOUS SPACES AND CLEAN ROOMS

Introducing the G2 VersaFlow Pump. A cutting-edge solution by G&G Technologies, this modular pump redefines the capabilities of pumps in Life-Sciences settings. Boasting a modular design, it accommodates an array of peristaltic and diaphragm pump heads for maximum adaptability. Engineered for HazLoc, ATEX, or general manufacturing space, Integrators will appreciate the versatility of the G2 VersaFlow Pump.

The G2 VersaFlow pump integrates with a variety of pressure and flow sensors for unmatched process control. An on-board data logging system provides valuable operational insights. G&G Technologies proudly presents the G2 VersaFlow pump as a symbol of reliability and innovation in fluid handling.



COMPLIES WITH:

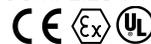

The First Modular SU Pump for Hazardous Locations

The G2 VersaFlow Pump addresses a critical need in the Single-Use Market: use in classified hazardous locations. Featuring an integrated purging system and rated components, the G2 VersaFlow pump is suitable for use in Division 2 and Zone 2 classified hazardous spaces. This approach maintains an IP 54 rating while pressurized, meeting standards and expectations for the life science industry.

Integration with a Broad Array of Sensors

The VersaFlow Pump may be integrated with an array of sensors to enable precise control over process parameters. Sensor options are available for general manufacturing spaces as well as hazardous locations. Sensor integration enables the pump to control a range of process parameters for a host of applications such as Nanofiltration, Depth Filtration, and Tangential Flow Filtration.




VF Pump Modules

The G2 VersaFlow Pump takes a universal approach to fluid transfer. A full array of popular pumps are configured into interchangeable modules which can be swapped in the field with minimal tooling. This industry leading versatility maximizes return on investment by eliminating the need for multiple pump systems from different manufacturers in one package. The G2 VersaFlow Pump offers a performance range from zero-flow up to 20 liters per minute. It's as easy as swap and go!

Fabricated from corrosion resistant materials including 316 Stainless Steel and polymers such as High Density Polyethylene and Acetal, the G2 VersaFlow Pump can withstand most industry cleaning agents



PUMP FLOW RATES L/min [L/hr]

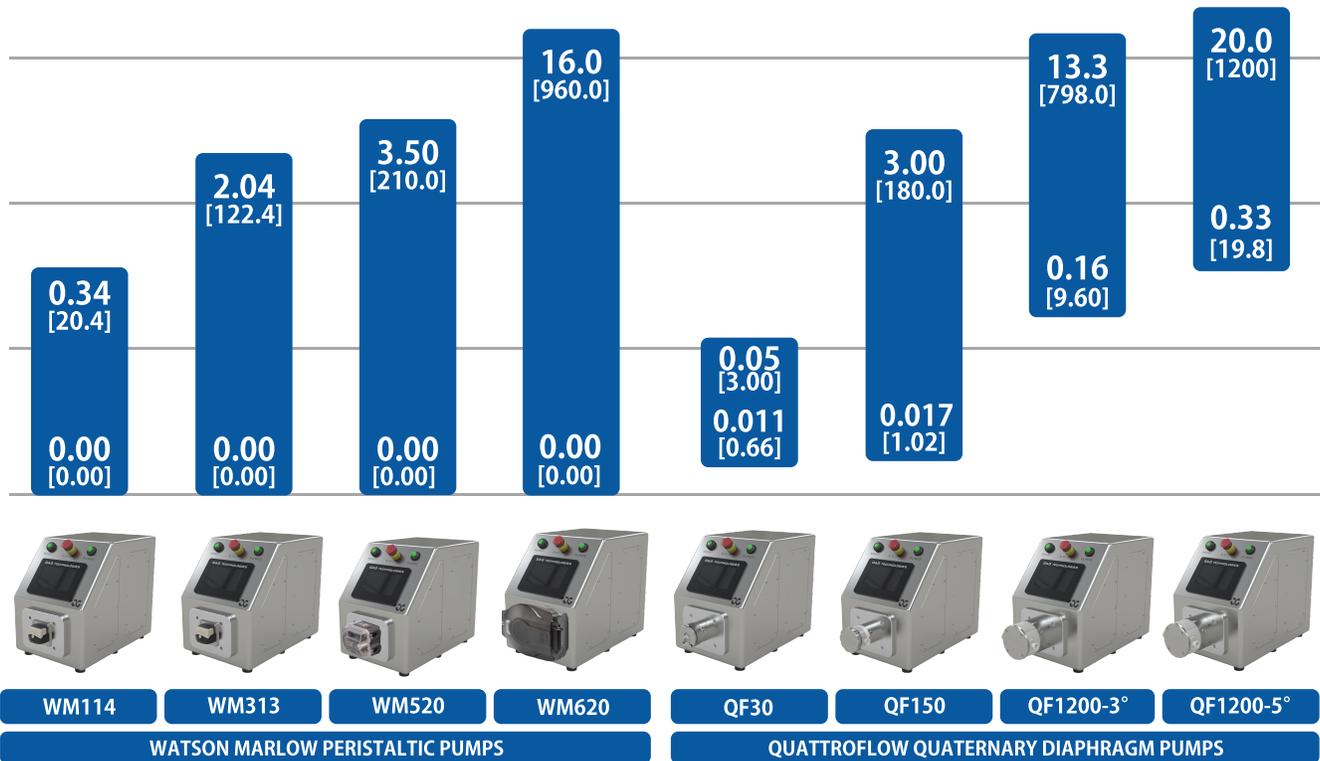

Contact G&G Technologies for customized solutions to your unique fluid transfer challenges

