



WORKS FASTER.

LASTS LONGER.

Medline OPA/28
High-level disinfectant solution.

SAFETY COUNTS AND MINUTES MATTER.

Time and safety are precious in processing. Getting more done in less time and reducing disinfectant exposure just got easier with Medline OPA/28. It works faster and lasts longer than most OPA brands to provide safety, convenience and value.

WORKS FAST.

Destroy pathogens in just 10 minutes (at room temperature) with this powerful, ready-to-use disinfectant. This reprocessor of semi-critical medical devices uses the active ingredient, 0.575% ortho-phthalaldehyde (OPA) to safely disable harmful pathogens, like:¹

- » MRSA
- » VRE
- » TB
- » HIV
- » Hepatitis Viruses




Bactericidal, fungicidal and tuberculocidal action in 10 minutes².



INCREASES CLEANING CYCLES.

Regain an extra disinfection cycle per hour with Medline OPA/28's fast cycling time². This can expand your operational efficiencies and reduce medical instrument downtime.

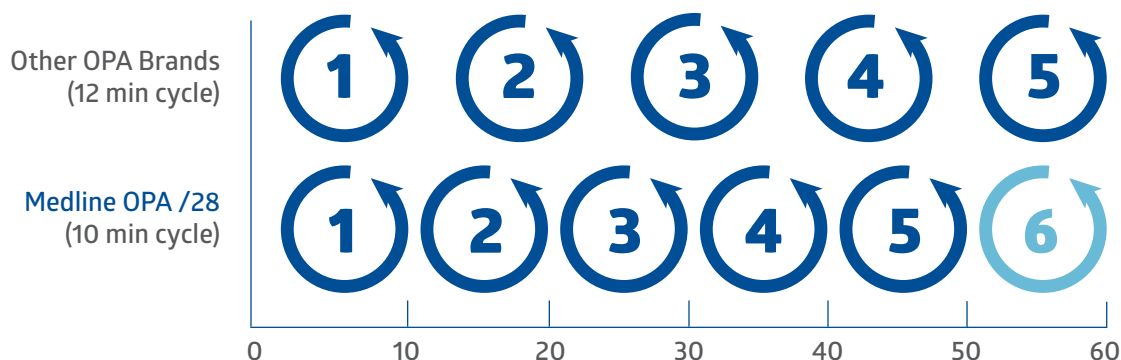

LASTS LONGER.

Medline OPA/28 has a reuse period of 28 days—up to twice as long as other OPA brands. One bottle can do the work of two over a 28-day period, creating cost-savings and safely limiting worker exposure to disinfectants.²



MEDLINE OPA/28 BY THE NUMBERS.

5 MINUTES

AER
DISINFECTION
TIME AT 25°C

10 MINUTES

MANUAL
DISINFECTION TIME
AT 20°C

28 DAY
REUSE
PERIOD



2 FINAL
RINSES
NEEDED IN
MEDIVATORS®
AER

3 FINAL
RINSES
NEEDED
FOR MANUAL
APPLICATIONS

GET FAST MRC VERIFICATION.

Test Strips are designed specifically for Medline OPA/28.

- » **Fast.** 3-second dip time; 90-second read time
- » **Intuitive.** Color coding provides pass/fail reference
- » **User friendly.** Large, easy-to-read indicator pad



References:

1. Data available upon request. 2. Provided that MRC meets requirements and product is used for full 28-day period.

MEDLINE OPA/28 HIGH-LEVEL DISINFECTANT

Discover an excellent way to reprocess heat-sensitive, semi-critical medical devices that are unsuitable for sterilization.

Item Number	Description	Pkg.
MDS88OPA28	High-level Disinfectant, 1-gallon bottles	4/cs
MDS88OPASTP	Test Strips, 50 strips/bottle	2/cs



Learn more about Medline OPA/28 by contacting your Medline representative, call 1-800-MEDLINE or visit us at www.medline.com



Medline Industries, Inc. One Medline Place, Mundelein, IL 60060

Medline United States
1-800-MEDLINE (633-5463)
medline.com | info@medline.com

Medline Canada
1-800-396-6996
medline.ca | canada@medline.com

Medline México
01-800-831-0898
medlinemexico.com | mexico@medline.com

FOLLOW US [f](#) [in](#) [t](#) [BLOG](#)