

CRONUS 8



PRINCIPLE

Cronus8, the latest innovation in skin rejuvenation technology, combines fractional microneedling with radiofrequency (RF) energy. Cronus8 offers patients a highly effective, minimally invasive solution for comprehensive skin rejuvenation.

Cronus8 delivers superior outcomes with a technology that aligns with the highest standards of care as to stay at the forefront of aesthetic



KEY BENEFITS FOR YOUR PRACTICE

Advanced Fractional Technology

Cronus8 penetrates deeper than traditional microneedling, delivering RF energy to the subdermal layers. This stimulates collagen production and remodels tissue for noticeable skin tightening and rejuvenation.

Versatile Treatment Areas

Cronus8's ability to address a wide range of skin concerns across multiple areas, including the face, neck, and body.

Minimal Downtime

A treatment that fits seamlessly into patients' busy lives. Cronus8 requires minimal recovery time, allowing an easy return to daily activities swiftly.

Optimized Patient Satisfaction

The precision and depth control of Cronus8 ensure tailored treatments that meet the unique needs of each patient.

KEY ADVANTAGES OF CRONUS 8

Dual-Layer Treatment

Effectively treats both the epidermis and dermis, ensuring comprehensive skin rejuvenation.

Versatile for All Skin Types

Safe and effective regardless of skin type, making it a universal solution.

Minimal Downtime

Patients experience quick recovery times, allowing them to return to daily activities sooner.

Innovative Database Technology

Streamlines treatment processes and enhances patient outcomes.

Reduced Risk of PIH

Lower risk of post-inflammatory hyperpigmentation, ensuring safer treatments.

Dual Fractional RF Modes

Offers both monopolar and bipolar RF modes for tailored treatments.

Accelerated Healing

Faster healing times with less discomfort, enhancing patient satisfaction.

PRECISION AND SAFETY

Adjustable Needle Depth

The needle depth (0.5-3.0mm) is precisely controlled by the system, protecting the handles from damage and ensuring consistent treatment results.

Stepping Motor Technology

Cronus8 utilizes a stepping motor to drive the needles accurately, ensuring controlled and stable movements. This reduces the risk of bleeding and enhances patient comfort.

WHY CHOOSE CRONUS8?

Non-Insulated vs. Insulated Tips

Non-Insulated Tips

Energy is distributed along the entire length of the needle, affecting both the dermis and superficial skin. This can result in longer recovery times due to the larger treated area.

Insulated Tips

Only the top 0.2mm of the needle delivers energy, minimizing damage to the superficial skin. This results in smaller wounds, less redness, and quicker recovery. The gold-coated stainless steel needles are highly durable, with a minimal diameter of 0.17mm, reducing pain and speeding up recovery.



VERSATILE TREATMENT OPTIONS

Acne Scars Treatment

Cronus8 effectively treats acne scars by targeting the dermal inflammatory response, stimulating healing, and reducing the appearance of scars. Bipolar RF technology ensures focused treatment with minimal discomfort.



Scar Reduction

RF energy penetrates deep into the skin, remodeling collagen to improve the appearance of scars over time. The result is smoother, more even skin.



Fine Wrinkles

Non-invasive RF treatments with Cronus8 tighten skin, stimulate collagen production, and rejuvenate the complexion without damaging the epidermal layers.




Skin Rejuvenation

Cronus8 promotes collagen, elastin, and hyaluronic acid production, reducing fine lines and wrinkles for a more youthful appearance.



Stretch Marks

RF skin resurfacing with Cronus8 delivers heat energy to stimulate collagen production, effectively treating and reducing the appearance of stretch marks on various body areas.

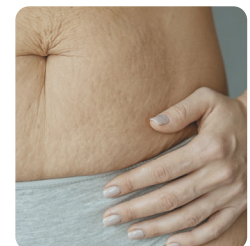


Pore Reduction

RF Microneedling with Cronus8 tightens skin and reduces pore size, providing a safe, year-round treatment option that's effective for all skin types.




HANDPIECES

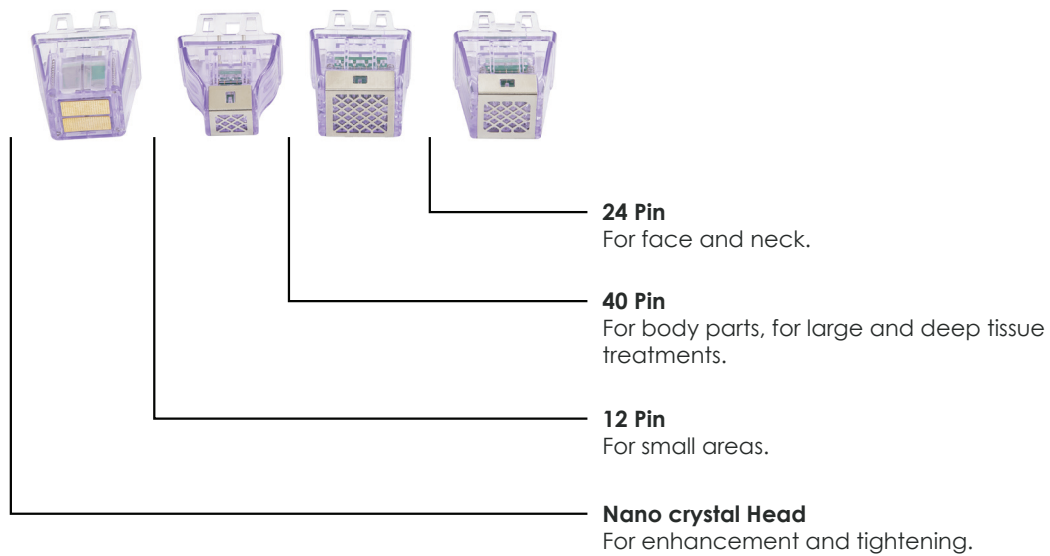

SPECIFICATIONS

TECHNOLOGY	Microneedle & RF
Touch Screen	15 inches
Therapy Probe	12p, 24p, 40p, Nano crystal, Thermal RF Heads
Depth	0.5 - 7mm
Operating Voltage	AC 110-220V \pm 10%, 50Hz - 60Hz
PACKAGE SIZE	103*50*50cm
GROSS WEIGHT	28Kg
NET WEIGHT	13.5Kg

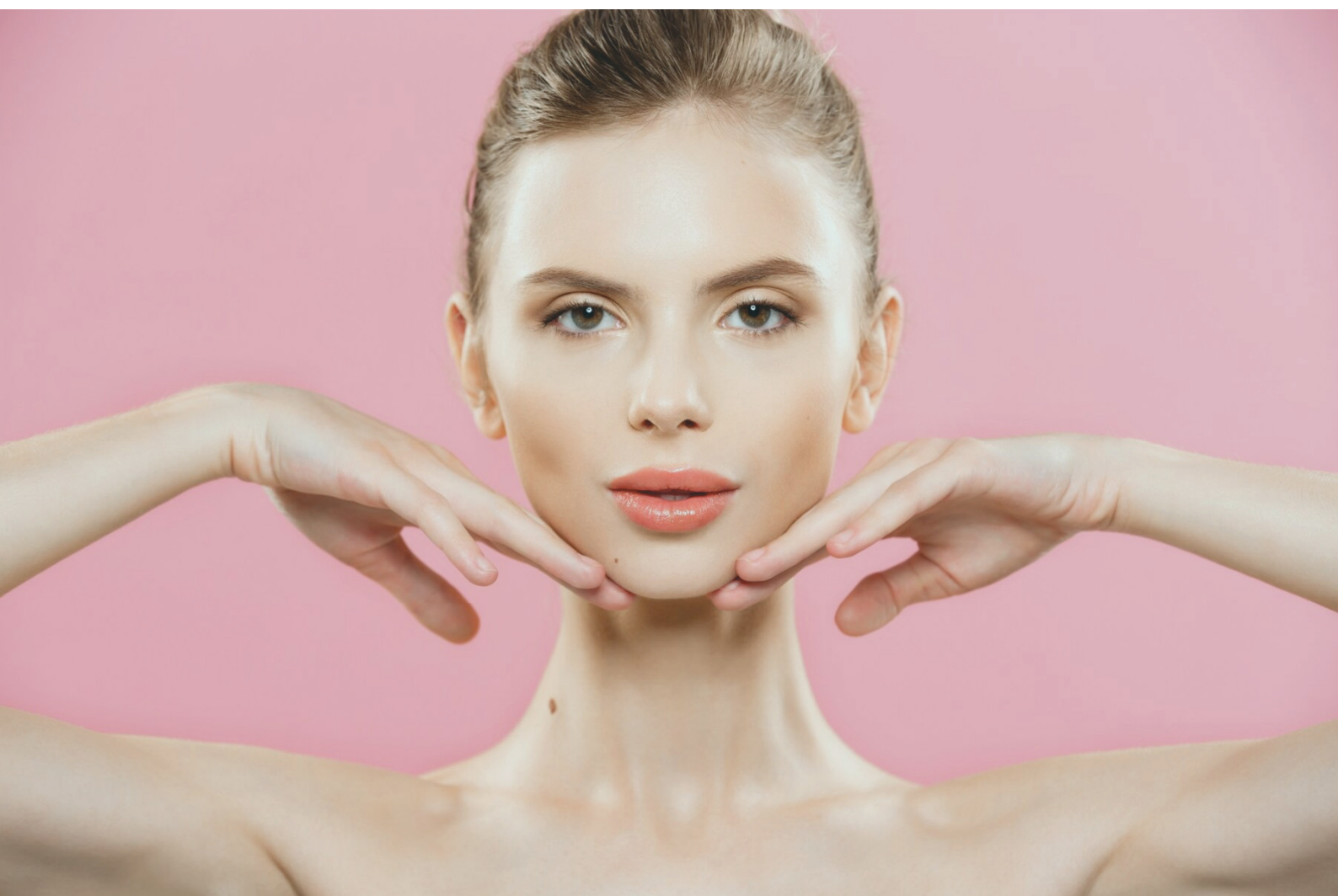



**Tomorrow's Technologies.
Today.**

www.axicontechnologies.com

© 2024 Axicon Technologies LLC - USA
All rights reserved.